

Mini-Bag Plus



1. Using your fingers, break the light blue seal between bag and medicine vial.



2. Holding the bag up, squeeze and release the bag multiple times to push fluid into the vial. Fill vial $\frac{1}{2}$ way, mix by gently swirling vial until all powder is liquid. Do not shake.



3. Holding medicine vial up with bag below, squeeze and release the bag, pushing air into the vial, until vial is empty. Check to make sure all powder is mixed and out. If residue is left, repeat steps 2 and 3.



4. Pull dark blue "stretchy" plug out of bottom of the mini-bag.

Mini-Bag Plus


5. Take tubing and clamp, then remove spike cover, push spike into bag. Twist back and forth, so it goes all the way in.



6. Hang bag on IV pole, squeeze chamber 1-3 times, until about 1/3 full. Resume with priming the tubing (see attached sheet for priming the tubing).




7. Set dial on tubing to desired rate. Refer to medication label/physician order for rate.

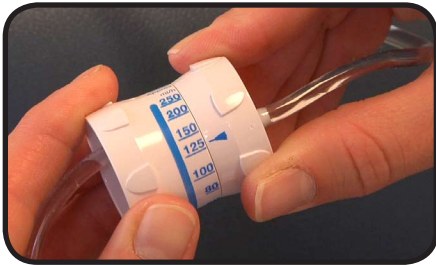

8. Once infusion is complete, make sure to keep fluid in the chamber when bag is empty. Clamp tubing (white clamp), disconnect tubing, and place a red cap on the end of tubing to keep the end from getting contaminated. Keep tubing on IV pole for the next dose (if within 24 hours), otherwise throw away.

Tubing needs to be changed every 24 hours!

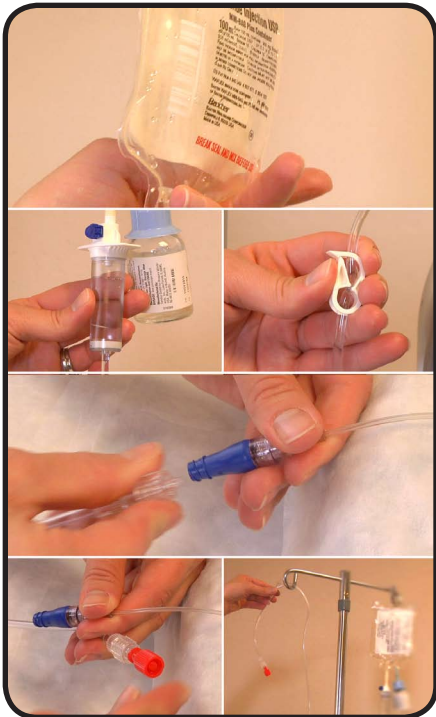

For any questions please contact Altius Healthcare

Altius
Healthcare
Infusion Therapy & Specialty Pharmacy
1-800-531-2469