



Priming Tubing

1. Open package and close white clamp and open the rate control wheel/dial to "open". When opening the package, if the dial is on "open," turn dial to close, then back to open to break the seal.

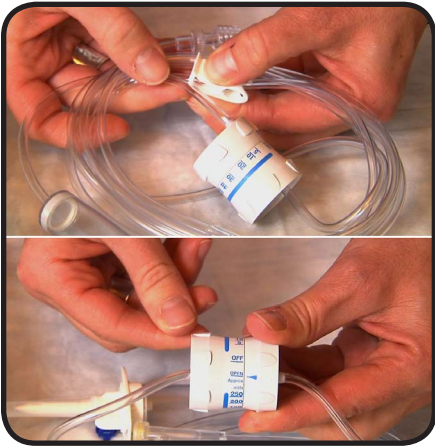


2. Remove the protective cover from tubing set's spike and medication bag, be careful not to touch spike to skin or other surfaces (keep sterile).

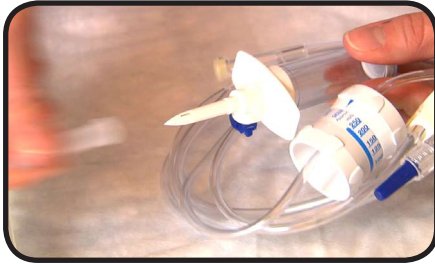

3. Insert the spike on the tubing set into the medication bag/container.



4. Hang container on the IV pole.

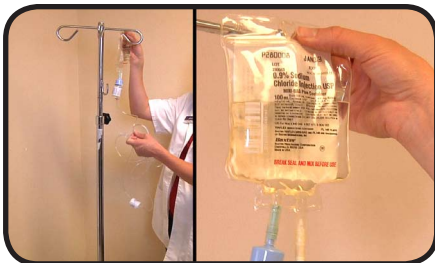


Priming Tubing


5. Prime the tubing by first squeezing the drip chamber to fill about 1/3 full.



6. Slowly open the clamp; hold the end of the tubing upright.

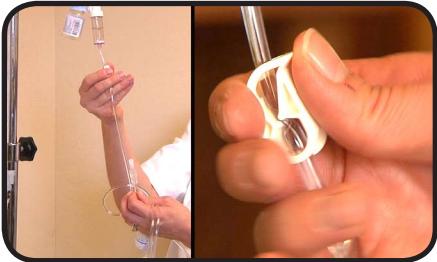

7. Prime entire length of tubing and then clamp. Look to make sure all the air is out of tubing (small intermittent bubbles are ok). If large air bubbles are in the tubing, call Altius Healthcare for assistance.

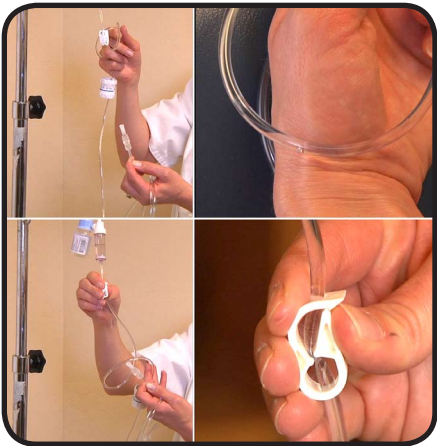

8. Remove protective cover from the end of the tubing set and connect to your IV line.



9. Set the rate control dial to the prescribed rate and unclamp tubing.

For any questions please contact Altius Healthcare

Altius
Healthcare
Infusion Therapy & Specialty Pharmacy
1-800-531-2469